



freedom shield

Direct Splash Face Protection

Neo One Distribution Ltd.

Unit 4c, Sett End Road North,
Shadsworth Business Park,
Blackburn, Lancashire,
United Kingdom, BB1 2PT

t. 01254 299888
e. sales@neoone.co
www.neoone.co

Company Registration No. 06882491
Company VAT No. GB 272805301


With a supply capacity of over 100,000 units per week. NeoOne are pleased to introduce the "freedom shield". This UK manufactured face shield is recommended for face protection against liquid droplets and splashes. The shield is not impact resistant and should not be used as protection against open flame or sparks.

The foam head spacer provides comfort for the user and a space for users to wear glasses if required. The adjustable head strap fits all head sizes easily and quickly. Includes User ID and Used Date entry location, for use monitoring and protection.



CE Certified

- Shield Made In The UK
- In Stock
- Supply Capacity in excess of 100,000 per week !
- Regulation (EU) 2016/425 Tests Passed
BS EN 166:2002 – 6.1, 6.3, 7.1, 7.1.1
EN 167:2001 – 3.2, 5
EN 168:2001 – 10.2, 10.2.2.2, 17, 19
Test report on reverse.

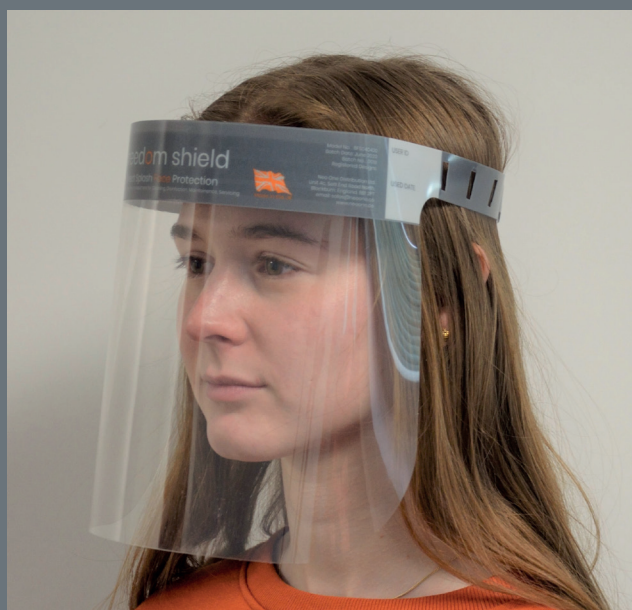

Pricing: June 2020

1 Box of 70no. shields
1 Pallet of 1,680no. shields
30 Pallets of 50,400no. shields
60 Pallets of 100,800no. shields

Carriage Inc. to a single UK mainland Address

£274 + VAT (price per shield £3.92)
£4,838+ VAT (price per shield £2.88)
£129,528 +VAT (price per shield £2.57)
£252,000 +VAT (price per shield £2.50)



freedom shield

Direct Splash Face Protection

Part No. : **NBFS040420**

Pallet Size: 1200mm wide x 800mm deep x 1118mm high (Euro Pallet)

Pallet Weight: 105kg. Number of face shields per Pallet: 1680no. Packed in 24no. boxes of 70 face shields.

Individual Box Weight 3.3kg

Individual Box Size: 380mm wide x 290mm deep x 305mm high.

Instructions: Store in a hygienic, dry, flat place.

Peel off selfadhesive backing & stick foam to rear of shield inside corner marks.

Insert tail through sizing hole. Hook head strap at preferred size to fit your head above your ears without contact with your ears for comfort.



Certificate Number: CE729963

PPE for use by healthcare professionals as per Commission recommendation 2020/403

Test Requirements.

BS EN 166:2002 + Results Table – Personal Eye Protection



Clause	Requirments	Result
6	Design and Manufacturing Requirements	
6.1	General construction Eye-protectors shall be free from projections, sharp edges or other defects which are likely to cause discomfort or injury during use.	Pass
6.3	Headbands Headbands, when used as the principal means of retention, shall be at least 10mm wide over any portion which may come into contact with the wearer's head. Headbands shall be adjustable or self-adjusting.	Pass
7	Basic, particular and optional requirements	
7.1	Basic requirements	
7.1.1	Field of vision The size of the field of vision is defined in conjunction with the appropriate head-form described in clause 17 of EN 168:2001.	Pass
7.1.2.1.2	Mounted oculars and unmounted oculars covering both eyes The refractive power characteristics of mounted oculars or unmounted oculars covering both eyes shall be measured by the method specified in clause 3.2 of EN 167:2001 at the visual centre of the ocular.	Pass
7.1.3	Quality of material and surface Testing shall be in accordance with clause 5 of EN 167:2001. Except for a marginal area 5 mm wide, oculars shall be free from any significant defects likely to impair vision in use, such as bubbles, scratches, inclusions, dull spots, pitting, mould marks, scouring, grains, pocking, scaling and undulation.	Pass
7.2.4	Protection against droplets and splashes of liquids b) face-shields cover the eye-region rectangle of the appropriate head-form as described in clause 10.2.2.2 of EN 168:2001 as assessed in accordance with 10.2 of EN 168:2001.	Pass
7.2.8	Lateral Protection Testing shall be in accordance with clause 19 of EN 168:2001.	Pass